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**Low/Medium Frequency  
Loin Massager**

- Product components:
1. Device
  2. Handheld Remote
  3. Power supply unit
  4. Electrode pads

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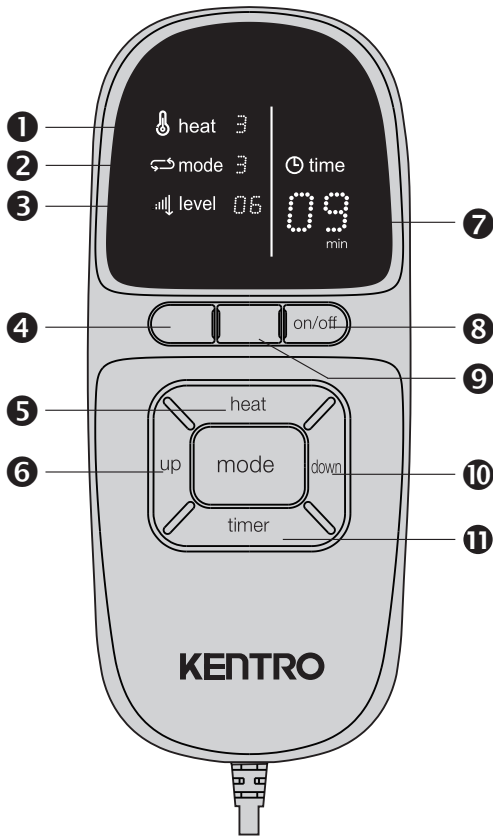
**Preparations**

1. Place the device on the bed, then connect the handheld remote and the mains plug to the device. After this, plug the mains plug into a socket.
2. Clean your waist with a moist cloth and lie down on the device with your back.

**Important! Both electrode surfaces have to be in direct contact with your skin. Otherwise you will not be able to operate the device.**

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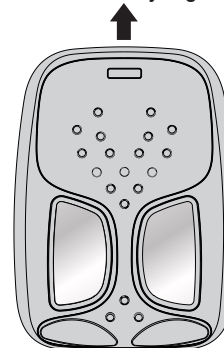
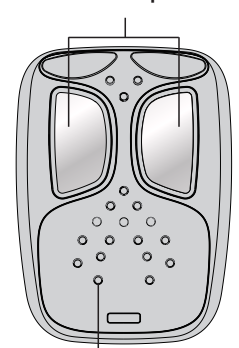
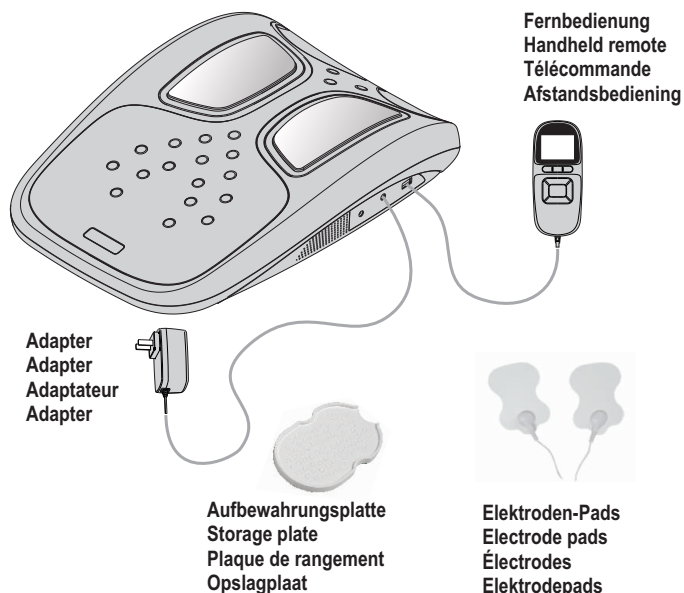
**Instructions for the use of the electrode pads**


1. Connect the electrode pads with the device via the cable.
2. Pull the protective strips off the electrode pads and stick the electrode pads where you need them.
3. After use, store the electrode pads on the storage plate in a dry, well ventilated place.
4. The surfaces of the electrode pads have to be clean. If they are not, please rinse them with clean water, then leave them to dry in the air. Each electrode can be used approx. 80 times (with an application of 20 minutes each). In case the sticking strength of the electrode diminishes even though they are cleaned, please replace them by new ones. This device only works with electrodes of the same type (55 x 80 x 2.5 mm). Please contact the specialist dealer you purchased the device from.

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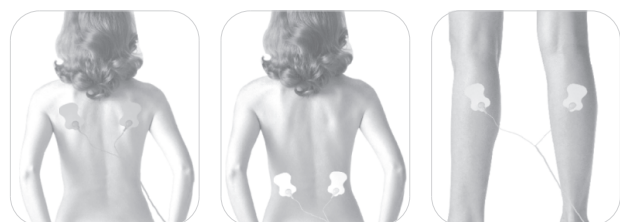The different modes simulate different massage techniques (patting, stroking, kneading). The first two modes provide a massage cycle of 3 different techniques.

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**Suggestion for the application:** Twice daily; massage function: 10-20 minutes; warming function: 15-30 minutes

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
**Safety Precautions**

**Warning!**  
In case of improper use of this product there is a risk of accidents, damage, or severe injury.

**Contraindications**  
Do not use this device with a pacemaker, heart-lung machine, ECG device etc. This could lead to dangerous situations.  
The device might not work flawlessly near shortwave or microwave therapy devices (within a 1-metre distance).  
The risk of heart flutter increases when the electrode surfaces are used near the chest.

**Warning!**  
This device may not be used by pregnant women and persons with sensitive skin, during menstruation, in case of heart disease, abnormal blood pressure, malignant tumors, bleeding problems in the brain, or acute disease without consulting a physician.  
1. Do not use this device in the vicinity of the heart, head, mouth, pubic area, or in case of skin diseases.  
2. Persons with cognitive inability or persons that are insensitive to heat may not use this device.  
3. Do not use the device in places with high temperatures, electromagnetic radiation, or in inflammable or moist places.  
4. Do not use the device when taking a bath or sleeping.  
5. Do not disassemble, repair, or extend the device to avoid a risk of malfunction or electric shock.

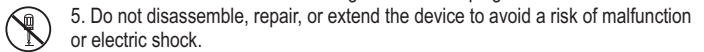**Caution!**  
1. Do not move the device during its operation. Switch it off before changing the area to be treated.  
2. This device may not be used on children or on persons who are not able to communicate.  
3. In case of complaints during the application of the device, immediately interrupt the application and see a physician.  
4. Unplug the device after each use or whenever the device is not in use.  
5. The device may not be used by persons (including children) with physical, sensory, or mental impairments or lack of experience or knowledge, unless they have been instructed to use the device or under supervision only.  
6. The device may be operated together with the accessories included in the delivery only.  
7. The device has a heatable surface. Persons insensitive to warmth should be careful when using the device.

For indoor use only.

Degree of protection: Type BF

Type of device: Type 2

The DC voltage supply anodes are inside the device and the cathodes on the outside.

**Cleaning and care**

- Cleaning with a dry, soft cloth.  
- After each use please clean the electrode pads with a moist, soft cloth.  
In case of strong soiling use a soft cloth moistened with medical alcohol (with a 75% concentration).  
Storage:  
- Keep out of reach of children.  
- Avoid places exposed to direct sunlight, high temperatures, or moisture.  
- Store in a dry, well ventilated place.  
Please keep the device out of reach of children in a dry and well ventilated place. The device should not be exposed to direct sunlight, strong heat, or moisture.

**Technical parameters and specifications**

1. Product specifications

Name: Low/Medium Frequency Loin Massager, model: 93204-KTR-301  
Power supply: 100-240 V AC; frequency: 50/60 Hz 0,5 A, which can be converted into 12 V DC, 1000 mA;  
equipment – safety category: Type BF application in category II  
Packaging dimensions 460 x 350 x 115 mm; weight: approx. 2450 g

2. Normal operating conditions

Environmental conditions: + 5°C to + 40°C; relative humidity: ≤80%; atmospheric pressure: 86-106 kPa

3. Electric power


Impulse frequency: 1-1200 Hz; pulse width: 20-400 µs; DC component: 0 V  
Electric energy of the single impulse with maximum output amplitude: >7 µC  
Maximum output power of the single impulse: ≤300 mJ  
Effective value of the maximum output amplitude: ≤25 V (50 mA)  
Output voltage peak during the measurement of the resistance (open circuit voltage): ≤500 V  
Output amplitude setting: steady and even with the minimum power of not more than 2% of the maximum output  
Load impedance: 500 Ω, with a ±10% tolerance; surface magnetic field strength: not less than 0,01 T  
Time setting: adjustable between 15 and 30 minutes, with a ±10% tolerance  
Warmth/temperature: +38°C + 48°C (with a temperature of 23°C ± 2°C)  
4. Transport and storage conditions  
Ambient temperature: -20°C to +55°C  
Relative humidity: ≤93 %  
Atmospheric pressure: 50 to 106 kPa

Problem	Check / Countermeasures	Solution
Handheld remote display does not work.	Check whether the cable of the handheld remote is correctly connected to the device.	Connect the cable correctly.
	Check the power connection to the socket.	Check the power connection to the socket.
	Check the power connection to the device.	Connect the cable correctly.
Handheld remote display works correctly but the user feels no effect.	The intensity button is set to 0 when the system is started.	Increase the intensity.
	No skin contact or insufficient skin contact of the device.	Put off clothes to ensure skin contact.


Sharp pain and redness of the skin.	The skin is too dry.	Moisten the skin.
	The electrode pads are soiled.	Clean the electrode pads as indicated in the instructions for use.
	Too long duration of application.	Recommended length of treatment: 15-20 minutes.
	Too high intensity	Reduce intensity level.
	The skin is too sensitive and does not respond to contact with the electrode pads.	Have there been allergies in the past? In case of mild allergies shorten the length of treatment. In case of serious allergies immediately interrupt the treatment.
	Cold or dry outdoor temperature could be the cause.	Either slightly moisten the skin before the application or use the warmth function for 3-5 minutes.

Problem	Check / Countermeasures	Solution
Weak stimulation	Insufficient skin contact with the electrode pads	Ensure direct skin contact and adapt the intensity and the massage mode
	Soiled electrode pads	Clean the electrode pads as indicated in the instructions for use.
No sensation of warmth	Warmth button switched off when the system is started.	Switch on the warmth button
No stimulation by the electrode pads.	The cable of the electrode pads is not correctly connected to the device.	Connect the cable correctly.
	The button is set to 0 when the system is started.	Set the button to 1 or above.


### 1. Lying on your back
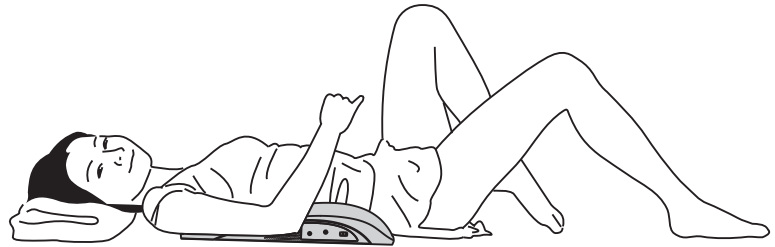
### 2. Getting up: After use, turn on your side and get up slowly with the help of your arms.
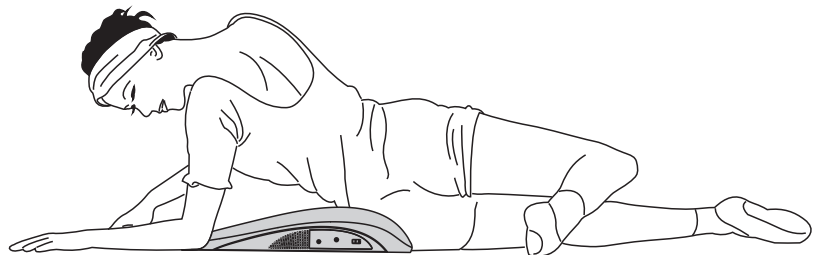
### 3. Do **NOT** raise your body abruptly.

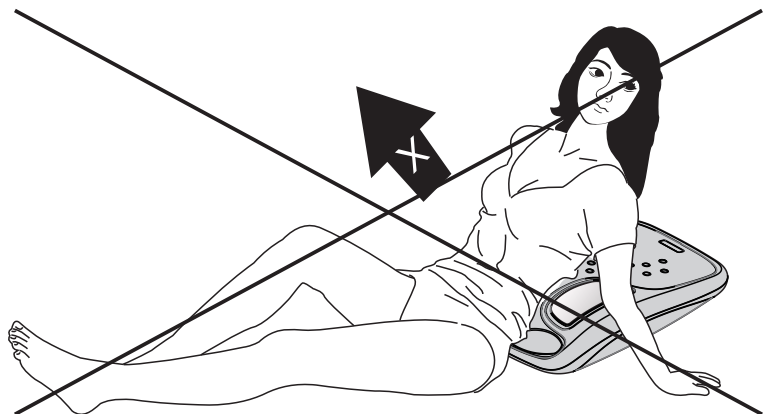**DISPOSAL:** Do not dispose of the device into ordinary household waste. Leave the device at a collection point for discarded electric devices instead.

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