GYMFORM® SIX PACK is a registered EU/CTM trade mark.  
International patents pending. All imitations will be prosecuted.

# GYM<sup>®</sup> form SIX PACK SPORT

## *Instruction Manual*

-  *Manual de instrucciones*
-  *Istruzioni per l'uso*
-  *Mode d'emploi*
-  *Manual de instruções*
-  *Bedienungsanleitung*
-  *Gebruiksaanwijzing*

Fig. 1 / Abb. 1

Fig. 2 / Abb. 2


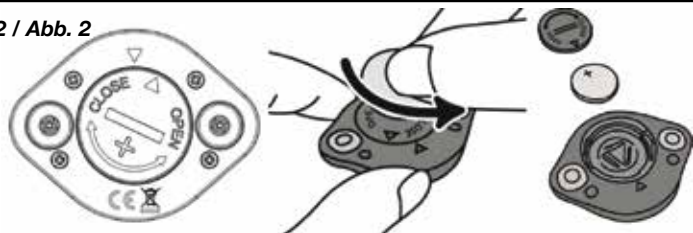
Fig. 3 / Abb. 3

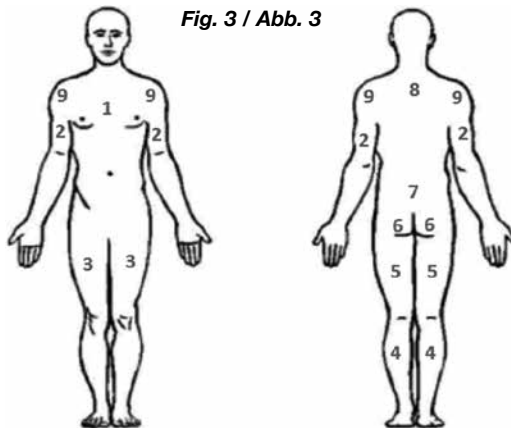Congratulations on your purchase of the Gymform® SIX PACK. Your safety is our primary concern, so please ensure that you familiarize yourself with all the information in this manual before using your Gymform® SIX PACK. It is essential that anyone who intends to use this device reads and understands the instructions provided. An informed user is a safe user. Thank you for purchasing the Gymform® SIX PACK. We hope you enjoy your sessions.

### WARNING

Before embarking on any exercise programme, including the Gymform® SIX PACK, consult your general practitioner and your physiotherapist.

Ensure that you read the instructions, warnings and precautions carefully before using your Gymform® SIX PACK. Gymform® SIX PACK is designed for healthy adults. However, the device should be used with care. Read the following guidelines and instruction manual before using the product.

- Applying electrodes in the thorax area may increase the risk of cardiac fibrillation.

- Using the Gymform® SIX PACK while connected to a high-frequency surgical unit may cause burns at the contact points of the Gymform® SIX PACK and may result in damage to your Gymform® SIX PACK.

- Operation of a short-wave therapeutic device or microwave oven in the general vicinity (e.g. 1 metre away) may cause instability in the output power of the Gymform® SIX PACK.

- This is a powerful electrical muscle stimulation device and may cause injury if you start with the most advanced exercise programme right away. Start with a lower level and gradually increase the intensity.

- Do not use on children or animals.

- Do not use while driving or sleeping.

- Do not use in the bathroom or in other places with a high amount of moisture or humidity.

- The batteries should be removed from the device before it is disposed of.

- The batteries should be safely disposed of.

- This device may be used by children aged 8 and over, by people with reduced physical, sensory or mental abilities, and by people lacking in experience and knowledge if they are appropriately supervised or given suitable instructions on the safe use of the device and understand the dangers its use entails. Children should not be allowed to play with the device. The cleaning and maintenance to be performed by the user should not be

done by children without supervision.

### DO NOT use Gymform® SIX PACK:

- If you have an electronic implant (such as a pacemaker) or suffer from another heart problem.

- If you are, or think you might be, pregnant.

- If you have cancer, epilepsy or are under medical supervision for a cognitive disorder.

- If you are breastfeeding.

- If you have a high fever.

### Women, before using Gymform® SIX PACK, you should wait:

• At least six months after the birth of your child (after consulting your doctor).

• A month after having an intrauterine device (IUD) fitted.

• At least three months after having a caesarean section (after consulting your doctor).

• Until after the first days of your menstrual period.

### Before using Gymform® SIX PACK you should obtain permission from your doctor or physiotherapist:

• If you suffer from a serious disease or injury not mentioned in this section.

• If you were operated on recently.

• If you are diabetic and take insulin.

### When putting on the product, always remember:

• Do not place the product on the neck, genitals, scars, moles, varicose veins or head.

• Do not place it on recent scars, inflamed skin, infected areas or areas that are susceptible to acne, thrombosis or other vascular problems (e.g. varicose veins), or any part of your body with limited sensitivity.

• Do not place on areas with injuries or limited mobility (e.g. fractures or sprains).

• Do not place the product directly on metal implants. Place receivers on the closest muscle.

### Possible Adverse Reactions:

There have been isolated reports of skin reactions caused by the use of electro-stimulation devices, including allergies, prolonged skin redness and acne. If any of these symptoms occur, switch off the device and consult your doctor.

**If you have never exercised before or have not exercised for a long time:**

Consult your general practitioner.

### Precautions:

• Do not use Gymform® SIX PACK if you are drowsy or while sleeping.

• Do not use Gymform® SIX PACK if your body is in contact with water.

• Keep the device out of reach of children.

• Do not use Gymform® SIX PACK in combination with any other device that transmits electrical current to your body or with any other machine.

• Do not connect the metal snap fasteners or electrodes to any other object.

• Do not modify the device.

• Use only the electrodes supplied by Industex.

• Do not attach two electrodes together.

• If any of the electrodes is not attached properly, you will notice an itching sensation.

• Do not use gel pads that are damaged, impaired or dirty, since these can cause burns or injuries.

• While using the product, make sure that the gel pads are in contact with the electrodes. If the electrode touches bare skin, it may cause injuries or burns.

• Attach the product to the area concerned and then switch the device on.

• When you have finished, switch the device off before detaching it from your body.

• If you experience any abnormal sensations on your skin or body, stop using it immediately and consult a doctor. Otherwise, you may have an accident or suffer an injury.

• If you use the product on full power or for a long period of time, it may cause muscle strain and tingling or redness of the skin.

• Remove jewellery and accessories such as necklaces, belts and earrings before using the device. Otherwise, you may have an accident or suffer an injury.

• Do not use the product with electronic devices such as watches, pedometers, mobile phones, etc., since these may affect its operation.

• Do not use products such as body creams or lotions with the device, since these can reduce its adhesive properties.

• The gel pad adheres to the skin better if hair has been removed from the area. However, do not use it on the same day you have removed the hair.

• Do not heat or cool the product before use, since it could become deformed. Do not subject the device to sudden temperature changes while in use, for example by pouring hot water on it or using it in a sauna, etc.

• Clean the area where you want to use it.

• This product is not waterproof. Do not use it in locations with high humidity levels or where it could come into contact with water.

**The set includes:**

- 1 Gymform® SIX PACK console
- 1 Gymform® SIX PACK for the abdominal area
- 1 Gymform® SIX PACK for other areas (two electrodes)
- 1 Gymform® SIX PACK belt
- 6 Self-adhesive pads
- 1 Instruction manual
- 1 Diet booklet

**PRODUCT DESCRIPTION (Fig. 1):**

1. 1 Gymform® SIX PACK console
2. 1 Gymform® SIX PACK for the abdominal area
3. 1 Gymform® SIX PACK for other areas (two electrodes)
4. 1 Gymform® SIX PACK belt
5. Self-adhesive pads for the Gymform® SIX PACK
6. Self-adhesive pads for the two-electrode silicone unit
7. "+ / I" button: increase intensity or switch on
8. "- / O" button: decrease intensity or switch off
9. Battery cover
10. Pilot light
11. Quick opening and closing system

**Installing the battery:**

The Gymform® SIX PACK unit runs on one CR2032 battery (not included).

1. Turn the console (1) upside down to open the battery cover (9). Insert a coin in the slot on the battery cover (9) and turn anticlockwise (see Fig. 2).
2. Insert the battery so that the "+ / I" terminal faces upwards. Ensure it is firmly attached. (See Fig. 2)
3. Replace the battery cover (9) and use the coin to turn the cover (9) clockwise to close it (see Fig. 2).

**Using the Gymform® SIX PACK:**

1. Attach the console (1) to the Gymform® SIX PACK silicone unit (2) using the snap fasteners. Ensure it is firmly attached.
2. Remove the plastic backing that protects the self-adhesive pads (5).
3. Place the six self-adhesive pads (5) on the Gymform® SIX PACK silicone unit (2). Each pad should be positioned just above the curved grey line.
4. Remove all plastic backings from the other side of the self-adhesive pads (5) and attach the Gymform® SIX PACK silicone unit (2) to the abdominal area (on the stomach).
5. Press the "+ / I" button (7) once to switch the device on. The pilot light (10) will flash to indicate that the device is set

to intensity level 1.

6. Press the "+ / I" button (7) again to increase the intensity level up to a maximum of 10.
7. Press the "- / O" button (8) to decrease the intensity level. When the device is set to intensity level 1, pressing the "- / O" button switches the unit off.
8. We recommend that you start with low intensity levels and gradually increase the level as you become more familiar with the device.
9. To switch the unit off directly, press the "- / O" button (8) for two seconds.
10. Gymform® SIX PACK has only one automatic programme. The programme will switch off automatically after 20 minutes.
11. You can achieve significant results in just 10 minutes a day, but if you want even better results, use the 20-minute programme every day.
12. Gymform® SIX PACK will switch off automatically if the self-adhesive pads are not in contact with the skin.
13. When the battery needs changed, Gymform® SIX PACK will switch off automatically.

**Note 1:** Gymform® SIX PACK comes with larger, extra plastic backings to protect the self-adhesive pads once they have been placed on the unit. After using the product, place the large plastic backings over the self-adhesive pads (5) to extend their life. Do not remove the self-adhesive pads from the unit until they need to be replaced with new ones.

**Note 2:** To ensure that the product fits perfectly, you can use the Gymform® SIX PACK belt and fasten it above the Gymform® SIX PACK unit.

- 1- Put the Gymform® SIX PACK unit on and switch it on. Choose the desired intensity level.
- 2- Place the belt (4) around your waist so that the quick opening and closing system (11) is on the outside.
- 3- Secure the belt (4) using the quick opening and closing system (11). It should feel comfortable and the unit should be securely fastened to the abdomen.

**IMPORTANT:** Never remove the two round orange labels on the back of the Gymform® SIX PACK silicone unit. If any of the labels become detached from the metal pins, replace them. If either of the two metal pins is exposed, do not use the product until it is covered again.

**WARNING:**

• As the battery runs low, the electrical stimulation pulse will become weaker. If the battery runs out, the pilot light (10) will switch off and the buttons may stop working. Replace the battery.

**Using the two-electrode silicone unit:**

1. This unit has been designed so that the product can be used on muscles other than the abdominals.
2. To do this, place the Gymform® SIX PACK console (1) on the front of the two-electrode silicone unit (3) using the snap fasteners. Ensure that they are firmly attached.
3. Remove the protective plastic backings from one side of the two self-adhesive pads (6) and place them on the silicone unit. They should cover the grey area
4. Remove the protective plastic backing from the other side of the self-adhesive pads (6) and attach the product to the area you want to work (see Fig. 3). Replace the protective plastic backings when you have finished using the product.

1. Chest
2. Arms (biceps and triceps)
3. Upper thigh (front)
4. Calves
5. Upper thigh (back)
6. Glutes
7. Lumbar
8. Upper back
9. Shoulders

**WE RECOMMEND THAT YOU REPLACE THE ELECTRODES AFTER EVERY 30 USES. IF YOU NOTICE THE ELECTRODES HAVE PREMATURELY LOST THEIR ADHESIVE POWER, WET THEM WITH A DAMP SPONGE.**

**NOTE:** If the area where the pads are attached starts to feel warm, the pads have probably reached the end of their useful life. Switch off the Gymform® SIX PACK, take it off, fit new pads and then resume your session.

**TECHNICAL SPECIFICATIONS:**

Shape of initial curve	Biphasic square wave pulse
Pulse duration	320 µs
Pulse frequency	3.3 ~ 33.3 Hz
Output voltage	36 V±5 V Load resistance 1000 ohm
Input voltage	3 V
Duration of session	Automatic power off after 20 minutes
Intensity	Adjustable from 1 to 10
Mode	1 variable frequency mode
Operating conditions:	10°C~40°C (50°F~104°F); 30%~85% RH
Storage conditions:	-10°C~50°C (14°F~122°F); 10%~95% RH
Maximum power consumption	32mA

This symbol indicates that the device is protected against electrostatic discharge by BF-type double insulation.

ISO 7000-1641 User Guide

**CLEANING AND MAINTENANCE**

DO NOT WASH THE Gymform® SIX PACK UNIT: WIPE WITH A DRY OR DAMP CLOTH.

Remove the batteries from your Gymform® SIX PACK if you do not plan to use it for an extended period of time. MAINTENANCE AND STERILIZATION OF ADHESIVE PADS

- Sterilize by wiping with wet fingers.
- Do not immerse in water.
- Keep in a cool, dry place at room temperature.
- Do not remove the insulating label.
- After using the device, always replace the protective plastic backings on the adhesive electrodes.
- Store in the bag supplied.

**PRECAUTIONS WITH BATTERIES**

Ensure the battery is properly installed before turning the unit on.

Turn the unit off to replace the battery.

Remove the battery if you are not going to use the device for an extended period of time.

Take the necessary precautions when disposing of the battery. Do not expose it to fire.

**TROUBLESHOOTING**

Problem	Cause	Solution
THE PILOT LIGHT DOES NOT COME ON	<ul style="list-style-type: none"> <li>- The battery has not been installed.</li> <li>- The battery has not been installed correctly.</li> <li>- The battery has run out.</li> </ul>	<ul style="list-style-type: none"> <li>- Install the battery.</li> <li>- Check that the "+" terminal of the battery is facing upwards.</li> <li>- Replace the battery.</li> </ul>
THE PILOT LIGHT COMES ON BUT THE UNIT DOES NOT WORK	<ul style="list-style-type: none"> <li>- There is no contact with the skin.</li> <li>- The console is not properly attached to the Gymform® Six Pack unit.</li> <li>- The intensity level is set to 0.</li> <li>- The microprocessor has crashed.</li> <li>- The battery is running low.</li> </ul>	<ul style="list-style-type: none"> <li>- Check that the gel pad is in direct contact with the skin.</li> <li>- Check that the console is properly attached to the Gymform® Six Pack unit using the snap fasteners.</li> <li>- Press the "+ /  " button to increase the intensity.</li> <li>- Restart it by removing the battery and replacing it after one minute.</li> <li>- Replace the battery.</li> </ul>

THE BATTERIES DO NOT WORK PROPERLY	<ul style="list-style-type: none"> <li>- The microprocessor has crashed.</li> <li>- The battery is running low.</li> </ul>	<ul style="list-style-type: none"> <li>- Restart it by removing the battery and replacing it after one minute.</li> <li>- Replace the battery.</li> </ul>
------------------------------------	--	---

**DISPOSAL OF MATERIALS**

A symbol of a crossed-out wheeled bin means you should find out about and follow local regulations about disposing of this kind of product.

Do not dispose of this product as you would other household waste. Dispose of this device in accordance with the corresponding local regulations.

Electrical and electronic devices contain hazardous substances that can have harmful effects on the environment and/or human health and should be recycled properly. The battery should be disposed of separately. To protect the environment, all batteries, including rechargeable ones, should not be disposed of as you would other household waste. Take them to a disposal centre and ensure it complies with applicable legislation on battery disposal.

**Quality Warranty:**

This product is covered by a warranty against manufacturing defects subject to the time periods stipulated by the legislation in force in each country.

This warranty does not cover damages resulting from inadequate use, negligent commercial use, abnormal wear and tear, accidents or improper handling.

**Australia & New Zealand only**

This product has the benefit of certain Consumer Guarantees. These are prescribed by the Australian Consumer Law & the New Zealand Consumer Guarantees Act 1993 both of which provide protection for consumers. There is no express warranty for this product in Australia or New Zealand. The above paragraph refers to other countries.

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Industex GmbH
ISL ITALY S.r.l.
Best Direct B.V.
Best Direct International Ltd,